


MATTRESS APOLO

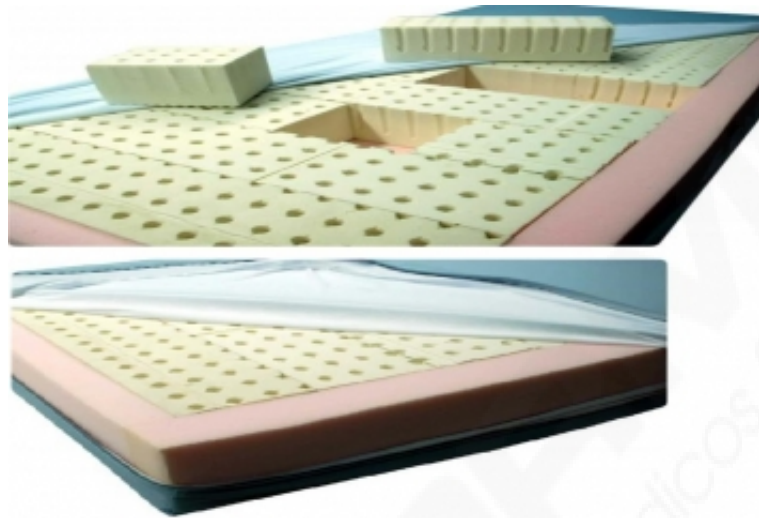

MAIN SPECIFICATIONS

- Mattress developed for the prevention and treatment of decubitus bedsores;
- It is a non-articulated and not reversible mattress;
- It is built with Latex cubes of 12 cm thick;
- Must remove one or more cubes Latex where it is intended to reduce the pressure made by the skin. In those places are placed lower cubes. These cubes come with this mattress. The gap has about 3 cm;
- The cubs network is supported by a base in Soyafoam with 5 cm thickness;
- Its constitution allows to minimize the pressure on the skin to avoid scarring and accelerate the healing of the existing bedsores;
- Ensures ergonomically correct body posture during sleep;
- Outer layer:
 - The mattress comes with a waterproof cover bielástica to liquids and gas-permeable, breathable that allows the breath of the interior, facilitating the hygiene of the mattress;
 - Fire-resistant, in accordance with BS 7175 1989, Crib 5;
 - The tissue can be washed to a temperature of 95°C with neutral detergent;
- Dimensions: 190x88x17cm;
- Note: we make mattresses on all measures (by budget).

Code	Description	Selling unit
61 CLAE AP	MATTRESS APOLO	CX. 1